

# Reliability and Validity of an Innovative Device for ACL Testing: The Mobil-Aider™

Dawn T. Gulick

**Background:** Knee disorders prevalence is estimated at more than 50% in a lifetime. There are over 250,000 anterior cruciate ligament (ACL) injuries per year in the United States. There are over 175,000 ACL reconstructions annually. This study was a double-blinded design to establish the reliability and validity of a new orthopedic device to measure linear translation of the tibia on the femur (ACL testing). **Methods:** A Zeiss Smartzoom microscope was used as the gold standard to assess the ability of the Mobil-Aider™ to measure linear translation. Sixty blinded measures were taken with each of 6 different devices. **Results:** Both the intraclass correlation and the Pearson correlation were .986. The Cronbach alpha reliability analysis was 0.992. Independent 1-sample *t* tests were performed on the differences between the Mobil-Aider™ and Zeiss values, and were not found to be significant ( $P = .42$ ); that is, they were the same. Bland–Altman plot and a linear regression revealed no propositional bias. Finally, with 360 measures over 6 devices, the power of this study was calculated to be 100%. **Discussion:** This data are the first step in establishing reliability and concurrent validity of a new device. As a result of the current data, the Mobil-Aider™ device is deemed a promising orthopedic tool for use in assessing the laxity of the ACL. Additional testing needs to be performed on both healthy and injured knees. **Conclusions:** There is potential for the Mobil-Aider™ to contribute to the assessment of ACL injuries, but additional human testing is needed.

**Keywords:** ACL injury, knee examination, knee stability testing

The prevalence of knee disorders is estimated at more than 50% in a lifetime.<sup>1,2</sup> There are approximately 250,000 to 300,000 anterior cruciate ligament (ACL) injuries per year in the United States.<sup>2,3</sup> There are over 175,000 reconstructions annually.<sup>2,3</sup> Athletes involved in high-demand sports have an increased risk of ACL injuries. Sports such as football, basketball, soccer, gymnastics, tennis, and skiing have been reported to elevate the risk of ACL injury.<sup>2</sup>

Although clinicians use a variety of examination techniques to assess the ACL,<sup>1</sup> there are 3 widely accepted clinical tests: anterior drawer, Lachman, and pivot shift.<sup>1,4–9</sup> These tests involve linear translation of the tibia on the femur. Meta-analysis and systematic reviews reported that these 3 tests have a wide range of diagnostic accuracy.<sup>10,11</sup> Individual modifications of the tests with varying degrees of tibial rotation and the magnitude of examiner experience have been reported to influence test accuracy.<sup>12,13</sup> The Lachman test is performed in the supine position with the knee in ~20° of flexion. Whether the test results are interpreted as dichotomous (positive = torn and negative = intact) or graded (grades I, II, and III), there is a subjective element to the data.<sup>13,14</sup> The clinician relies on the subjective assessment of end feel and being able to determine millimeters of translation to assign a grade of I (<5 mm), II (5–10 mm), or III (>10 mm).<sup>14</sup>


A study by Makhmalbaf et al<sup>5</sup> reported the sensitivity of the Lachman test to be 93.5% with a significant difference between men and women (66.7% vs 94.6%). A study by Guillodo et al<sup>15</sup> found the misdiagnosis of acute ACL injuries by emergency room physicians to be 74%. There are many reasons for these statistics. When there is a mismatch between the girth of the patient's leg and the size of the clinician's hand, it can be very challenging to

perform a Lachman test. If the patient's knee is in too much flexion or if the hamstring musculature is not relaxed, false-negative results may be obtained. Furthermore, asymmetry in side-to-side laxity or a soft end point is considered abnormal. Even for the most experienced clinician, a criterion of 3 to 5 mm is extremely challenging to quantify by "feel." Having an instrument to quantify the linear translation of this ligament would be very valuable.<sup>13,14</sup>

Magnetic resonance imaging (MRI) is often used for ACL injury diagnosis.<sup>3</sup> However, an MRI is a static image and is best used in conjunction with instrumented laxity devices to assess biomechanical behavior.<sup>16</sup> Furthermore, the ability of MRI to identify partial ACL tears has been called into question.<sup>17</sup> Over the past 3 decades, there have been a few instruments reported to be able to assess ACL laxity. However, all of the devices have been met with challenges. The Hall effect strain transducer was implantable.<sup>18</sup> The Rottometer was a computer-assisted goniometer used to measure rotation of the tibial axis.<sup>19</sup> The Vermont knee laxity device was very bulky and required a significant amount of time to utilize.<sup>20</sup> The Lars rotational laxiometer was dependent on too many variables.<sup>21</sup> The Kinematic Rapid Assessment requires Bluetooth technology to measure acceleration of the tibia on the femur.<sup>22</sup> The Telos is used in conjunction with radiographs.<sup>23–26</sup> The Vernier dial test indicator has not been validated.<sup>27</sup> The Telos, GNRB, and KT1000/2000 require considerable set-up time and do not involve direct clinician contact with the patient.<sup>28,29</sup> The GNRB sells for \$13,800, and the KT1000/2000 is no longer being produced (only available through the resale market). Several studies have reported substantial variability in the measures using the KT1000/2000, with a false-negative rate as high as 28%.<sup>27,30</sup> Despite these drawbacks, the KT1000/2000 device did sell because it met an unfulfilled need.<sup>27,30–32</sup>

To address this need, the Mobil-Aider™ (Therapeutic Articulations, LLC, Spring City, PA) was developed (Figure 1). The device is lightweight (<370 g; <13 oz) and used to measure linear

Gulick is with Physical Therapy, Widener University, Chester, PA. Gulick (dtgulick@widener.edu) is corresponding author.

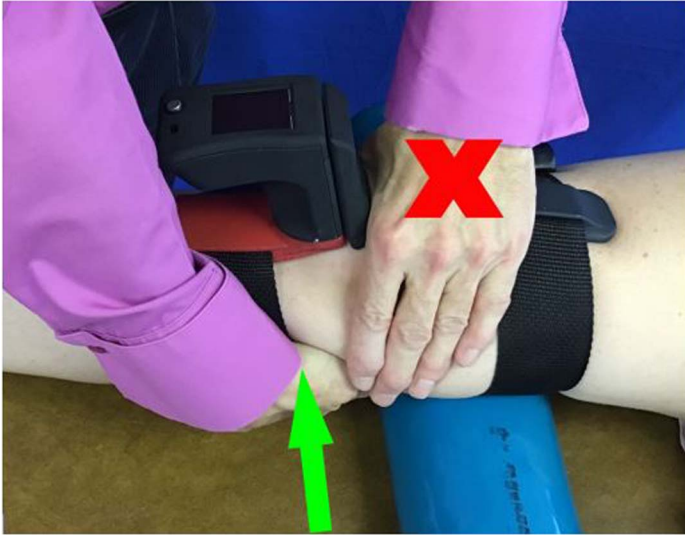

**Figure 1** — Lachman test with Mobil-Aider device.

translation of a joint. When the device is stabilized against the proximal joint surface with the axis aligned with the joint, the distal portion of the device is translated with the distal segment of the body. The device assesses the linear translation and displays it in a light emitting diode (LED) display in millimeters. The device permits the clinician to interface the device between the patient's skin and the clinician's hands, so there is no need to deviate from the standardized Lachman technique. In fact, the Mobil-Aider™ can even be used to perform a prone Lachman test. The contour of the device conveniently puts the knee in the optimal 20° to 30° knee-flexion position.<sup>9</sup> The nonelastic straps hold the device firmly in place to allow the clinician to attend to the technique and appreciate the qualitative end feel.<sup>33</sup> The LED reading (in millimeters) will allow for right and left comparisons, as well as serial measures of a given knee after surgery or throughout the rehabilitation process. The purpose of this study was designed to establish the reliability and validity of a new orthopedic device to measure linear translation of the tibia on the femur (ACL testing).

## Methods

The gold standard used to assess the Mobil-Aider™ was the Zeiss Smartzoom 5 Microscope (Carl Zeiss Microscopy GmbH, Jena, Germany). The Zeiss is designed for metallographic analysis; inspection of aerospace, automotive, and electronic components; and device failure analysis. Magnification ranges from 10× to 1011× with coaxial illumination. All data were collected in the University of Pennsylvania clean room. The Zeiss is self-calibrating. The metal measurement devices were positioned in parallel and secured to the Mobil-Aider™. The Mobil-Aider™ was positioned on the Zeiss stage. With the Mobil-Aider™ in focus, a screenshot of the baseline position of the device was obtained and saved as “Baseline1.” The Mobil-Aider™ device was translated to a random distance. The translation was revealed on the Mobil-Aider™ LED screen (in millimeters) and recorded on the data sheet. Another Zeiss screenshot was taken and saved as “M1.” The Mobil-Aider™ was reset, and the process was repeated for a total of 60× across 6 different serialized Mobil-Aider™ devices. Thus, a total of 360 measures were performed. The images saved from

the Zeiss were then reviewed, measured, and recorded on a data collection sheet. At the conclusion of all measurements, the data sheets were matched for Mobil-Aider™ and Zeiss measures. This resulted in all measures being blinded.

Intraclass correlation coefficient and Pearson correlations were performed to determine how strongly the measures of the 2 devices resemble each other. A Cronbach alpha reliability analysis was performed to measure internal consistency (a measure of how well a test addresses different constructs and delivers reliable scores). Independent 1-sample *t* tests were performed to determine if the 2 sets of data were significantly different from each other. A Bland–Altman plot was also generated, and a linear regression was calculated to check for propositional bias. Finally, a power calculation was performed to determine the confidence in the data.

## Results

The data analysis was performed with SPSS Statistics 23 software (IBM, Chicago, IL) and began with an intraclass correlation coefficient to measure the association of the measures as pooled means and SDs. Pearson correlation coefficients were calculated because it is deemed the best method of measuring the association between variables based on the method of covariance.<sup>34</sup> The Pearson correlation is a linear index. The correlations ranged from .986 to .997 and demonstrate a strong relationship between the 2 measures, but they do not confirm reliability or validity. Cronbach alpha is a measure of reliability. The analysis was performed on each device, as well as the overall measures of all devices. This was done because the alpha coefficient can be increased by simply increasing the number of items on the analysis. Cronbach alpha was ranged from 0.992 to 0.997. Independent 1-sample *t* tests were performed on the differences between the Mobil-Aider™ and the Zeiss values. This was performed to confirm that the measures were the same, that is, validity. The independent 1-sample *t* test for all devices was not found to be significant at the  $P = .05$  level ( $P = .42$ ). This indicates that the measures were not significantly different, that is, they were the same. In addition, the standard error of the mean was calculated because it is a measure of the dispersion of sample means around the population mean. A low value is a positive reflection of the accuracy of the data. A graph of the values is displayed in Figure 2. The Bland–Altman plot (Figure 3) displays the mean difference and the 95% upper and lower confidence limits. No propositional bias was identified. Finally, with 360 measures, the power of this study was calculated to be 100%.

## Discussion

The results of this study demonstrate both the reliability and the concurrent validity of the Mobil-Aider™ in a laboratory setting. This is the first step in assessing the clinical viability of any device. The ability to quantify ACL laxity could be a valid component to assessing the magnitude of an injury. Of course, the Mobil-Aider™ cannot identify the compromised bundles. However, with additional research, one may be able to relate the delta of the linear displacement of the involved versus the uninvolved knee to a grade of injury. Furthermore, the complete assessment of an ACL injury is based on the mechanism of injury, the clinical presentation, and imaging. The mechanism of injury can be contact or noncontact in nature. As the primary restraint to anterior translation of the knee,

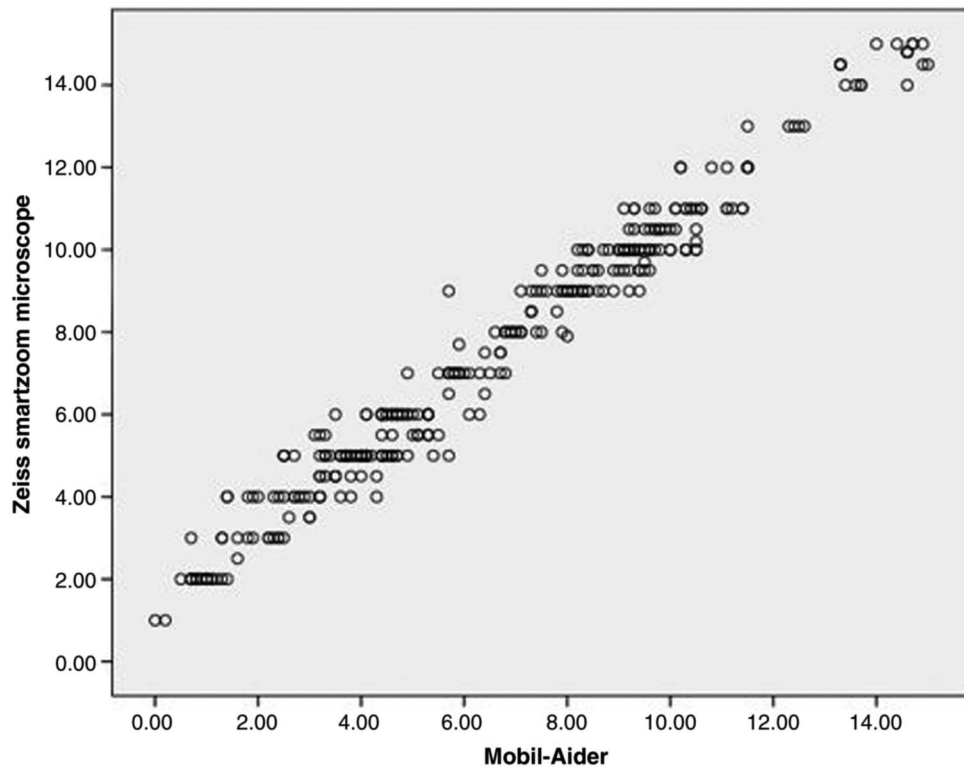


Figure 2 — Scatterplot of the Pearson correlation coefficient.

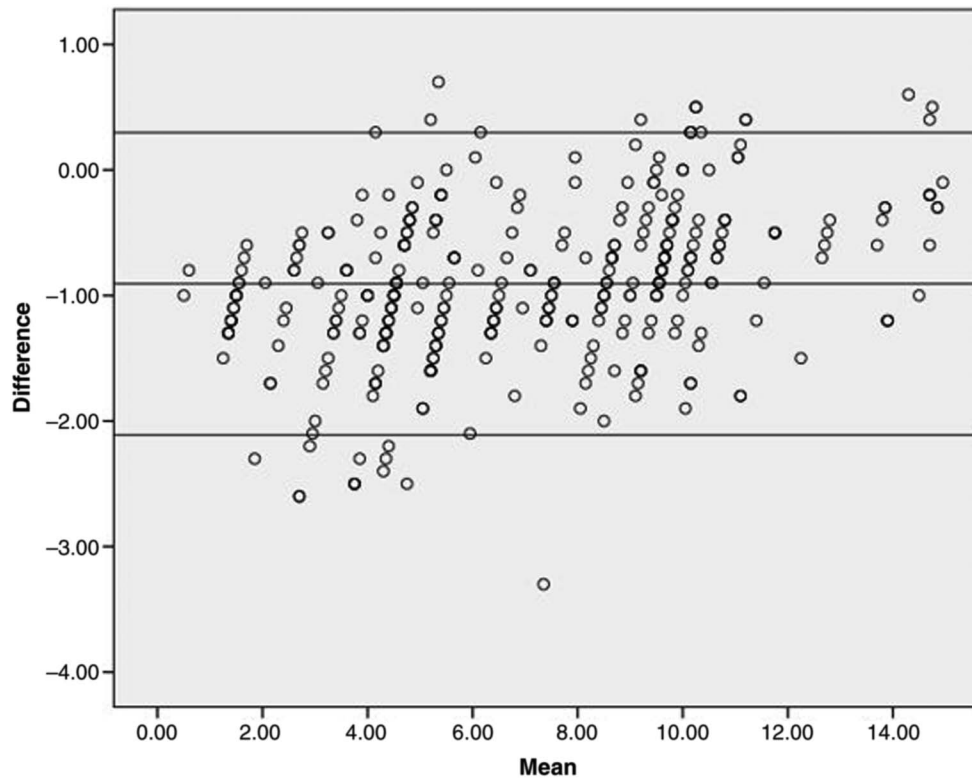


Figure 3 — Bland–Altman plot with 95% upper and lower confidence limits.

rapid deceleration, change in direction of movement, and knee hyperextension are just a few of the many ways the ACL can be damaged.<sup>14</sup> People report feeling a popping sensation and experiencing swift swelling with instability.

There are challenges in the clinical testing process: the subjective assessment of the Lachman test, the potential false negative if excessive knee flexion permits the posterior horn of the meniscus to influence anterior linear translation, and the mismatch of the size of the patient's limb to that of the clinician's hand can make translating the tibia on the femur very difficult. In addition, detecting millimeters of excessive translation while firmly grasping the limb can be difficult.

Magnetic resonance imaging is often used for ACL injury diagnosis.<sup>3</sup> These images only provide a "static" assessment of the ACL condition. It does not provide any information about the biomechanical behavior of the ligament or the joint.<sup>16</sup> Furthermore, the ability of MRI to identify partial ACL tears has been called into question by the research of Mikashima,<sup>17</sup> and the imaging process is very expensive.

The KT2000 device was available until 2012 but was large and cumbersome, with the need for 2 additional stabilizing components for proper positioning. The KT2000 used a handle to translate the tibia. The interface between the poorly contoured patella counterforce, the Velcro straps on the femur and tibia, and the handle through which the force was translated all contribute to potential error in the measurements. Several studies have reported substantial variability in the measures using the KT2000, with reports suggesting a false-negative rate as high as 28%.<sup>27,30–32</sup> Despite these drawbacks, the KT1000/2000 device met an unfulfilled need until 2012. Since that time, it is no longer on the market. The desire to address these clinical issues led to the development of the Mobil-Aider™ device to quantify tibial anterior translation with a simple handheld device that mimics the classic Lachman test. The contoured components allow the clinician to be in direct contact with the patient. The ability to provide a quantified, digital readout can provide significant assistance to the clinician assessing the ACL of the involved knee to that of the uninvolved knee.<sup>33</sup>

## Conclusions

The data collected in this study are the first step in establishing reliability and concurrent validity of a new device. As a result of the current data, the Mobil-Aider™ device is a promising orthopedic tool for use in measuring the linear translation of the tibia on the femur. Increased anterior translation of the knee may indicate underlying knee pathology, which may include an ACL injury. Additional testing will need to be performed on both healthy and injured knees. This is the next reasonable step.

## References

1. Decary S, Ouellet P, Vendittoli P-A, Desmeules F. Reliability of physical examination test for the diagnosis of knee disorders: evidence from a systematic review. *Man Ther*. 2016;26:172–182. PubMed ID: 27697691 doi:10.1016/j.math.2016.09.007
2. Texas Sports Medicine and Orthopedic Group. <https://www.txsportsmed.com/acl-tear-knee-specialists-dallas-frisco-fort-worth-tx/>. Accessed Nov 1, 2018.
3. Huang W, Zhang Y, Yao Z, Ma L. Clinical examination of anterior cruciate ligament rupture: a systematic review and meta-analysis. *Acta Orthop Traumatol*. 2016;50(1):22–31.

4. Galway RD, Beaupre A, MacIntosh DL. Pivot shift: a clinical sign of symptomatic anterior cruciate insufficiency. *J Bone Joint Surg*. 1976;54B(4):763–764.
5. Makhmalbaf H, Moradi A, Ganji SG, Omidi-Kashani F. Accuracy of Lachman and anterior drawer tests for anterior cruciate ligament injuries. *Arch Bone Joint Surg*. 2013;1(2):94–97. PubMed ID: 25207297
6. Malanga GA, Andrus S, Nadler SF, McLean J. Physical examination of the knee: a review of the original test description and scientific validity of common orthopedic tests. *Arch Phys Med Rehabil*. 2003;84:592–603. PubMed ID: 12690600 doi:10.1053/apmr.2003.50026
7. Marshall JL, Wang JB, Furman W, Girgis FG, Warren R. The anterior drawer sign: what is it? *J Sports Med*. 1975;3(4):152–158. PubMed ID: 1219194 doi:10.1177/036354657500300402
8. Peeler J, Leiter J, MacDonald P. Accuracy and reliability of anterior cruciate ligament clinical examination in a multidisciplinary sports medicine setting. *Clin J Sports Med*. 2010;20(2):80–85. doi:10.1097/JSM.0b013e3181ccca45
9. Torg JS, Conrad W, Kalen V. Clinical diagnosis of anterior cruciate ligament instability in the athlete. *Am J Sports Med*. 1976;4(2):84–93. PubMed ID: 961972 doi:10.1177/036354657600400206
10. Leblanc M-C, Kowalczyk M, Andruszkiewicz N, et al. Diagnostic accuracy of physical examination for anterior knee instability: a systematic review. *Knee Surg Sports Traumatol Arthrosc*. 2015;23:2805–2813. PubMed ID: 25763847 doi:10.1007/s00167-015-3563-2
11. Scholten RJ, Opstelten W, van der Plas CG, Bijl D, Deville WL, Bouter LM. Accuracy of physical diagnostic tests for assessing ruptures of the anterior cruciate ligament: a meta-analysis. *J Fam Pract*. 2003;52(9):689–694. PubMed ID: 12967539
12. Draper DO, Schulthies SS. Examiner proficiency in performing the anterior drawer and Lachman tests. *J Orthop Sports Phys Ther*. 1995;22(6):263–266. PubMed ID: 8580954 doi:10.2519/jospt.1995.22.6.263
13. Jarbo KA, Hartigan DE, Scott KL, Patel KA, Chhabra A. Accuracy of the Lever sign test in the diagnosis of anterior cruciate ligament injuries. *Orthop J Sports Med*. 2017;5(10):2325967117729809. doi:10.1177/2325967117729809
14. Koster CH, Harmsen AMK, Lichtenberg MC, Bloemers FW. ACL injury: how do the physical examination tests compare? *J Fam Pract*. 2018;67(3):130–134. PubMed ID: 29509817
15. Guillodo Y, Rannou N, Dubrana F, Lefevre C, Saraux A. Diagnosis of anterior cruciate ligament rupture in an emergency department. *J Trauma*. 2008;65(5):1078–1082. PubMed ID: 19001976 doi:10.1097/TA.0b013e3181469b7d
16. Espregueira-Mendes J, Andrade R, Leal A, et al. Global rotation has high sensitivity in ACL lesions within stress MRI. *Knee Surg Sports Traumatol Arthrosc*. 2017;25:2993–3003. PubMed ID: 27530386 doi:10.1007/s00167-016-4281-0
17. Mikashima Y. Can MRI distinguish between a partial anterior cruciate ligament (ACL) tear and a normal ACL? *Int J Orthop Rehabil*. 2016;3:51–57. doi:10.12974/2313-0954.2016.03.02.2
18. Howe JG, Wertheimer C, Johnson RJ, Nichols CE, Pope MH, Beynon B. Arthroscopic strain gauge measurement of the normal anterior cruciate ligament. *J Arthrosc Relat Surg*. 1990;6(3):198–204. doi:10.1016/0749-8063(90)90075-0
19. Almquist PO, Arnbjornsson A, Zatterstrom R, Ryd L, Ekdahl C, Fridén T. Evaluation of an external device measuring knee joint rotation. *J Orthop Res*. 2002;20:427–432. PubMed ID: 12038614 doi:10.1016/S0736-0266(01)00148-6
20. Uh BS, Beynon BD, Churchill DL, Haugh LD, Risberg MA, Fleming BC. A new device to measure knee laxity during



- weightbearing and non-weightbearing conditions. *J Orthop Res.* 2001;19:1185–1191. PubMed ID: [11781022](#) doi:[10.1016/S0736-0266\(01\)00055-9](#)
21. Bleday RM, Fanelli CG, Giannotti BF, Edson CJ, Barrett TA. Instrumented measurement of the posterolateral corner. *Arthroscopy.* 1985;14:489–494. doi:[10.1016/S0749-8063\(98\)70077-5](#)
  22. Berruto M, Uboldi F, Gala L, Marelli B, Albisetti W. Is triaxial accelerometer reliable in the evaluation and grading of knee pivot-shift phenomenon? *Knee Surg Sports Traumatol Arthrosc.* 2013; 21:981–985. doi:[10.1007/s00167-013-2436-9](#)
  23. Dejour DD, Ntagiopoulos PG, Saggin PR, Paneset J-C. The diagnostic value of clinical test, magnetic resonance imaging, and instrumented laxity in the differentiation of complete versus partial anterior cruciate ligament tears. *J Arthrosc Relat Res.* 2013;29(3): 491–499.
  24. Jardin C, Chantelot C, Migaud H, Gougeon F, Debroucker MJ, Duquennoy A. Low accuracy of the KT1000 arthrometer versus Telos radiographic measurements of anterior laxity after ACL graft. *Revue de Chirurgie Orthopedique.* 1999;85:698–707.
  25. Jorn LP, Friden T, Ryd L, Lindstrand A. Simultaneous measurements of sagittal knee laxity with an external device and radiostereometric analysis. *J Bone Joint Surg.* 1998;80:169–172. doi:[10.1302/0301-620X.80B1.7811](#)
  26. Lerat JL, Moyen B, Jenny JY, Perrier JP. A comparison of pre-operative evaluation of anterior knee laxity by dynamic x-rays and by the arthrometer KT100. *Knee Surg Sports Traumatol Arthrosc.* 1993;1:54–59. doi:[10.1007/BF01552160](#)
  27. Barcellona MG, Christopher T, Morrissey MC. Bench testing of a knee joint arthrometer. *Orthopedics.* 2013;36(8):e1000–e1006. PubMed ID: [23937744](#) doi:[10.3928/01477447-20130724-14](#)
  28. Ryu SM, Na HD, Shon OJ. Diagnostic tools for acute anterior cruciate ligament injury: GNRB, Lachman Test, and Telos. *Knee Surg Relat Res.* 2018;30(2):121–127. PubMed ID: [29554717](#) doi:[10.5792/ksrr.17.014](#)
  29. Vauhnik R, Morrissey MC, Perme MP, Sevsek F, Rugelj D. Inter-rater reliability of the GNRB knee arthrometer. *Knee.* 2014; 21:541–543. PubMed ID: [24239304](#) doi:[10.1016/j.knee.2013.10.004](#)
  30. Graham GP, Johnson S, Dent CM, Fairclough JA. Comparison of clinical tests and the KT100 in the diagnosis of anterior cruciate ligament rupture. *Br J Sports Med.* 1991;25(2):96–97. PubMed ID: [1751898](#) doi:[10.1136/bjism.25.2.96](#)
  31. Forster JW, Warren-Smith CD, Tew M. Is the KT1000 knee ligament arthrometer reliable? *J Bone Joint Surg.* 1989;71-B(5):843–847. doi:[10.1302/0301-620X.71B5.2584257](#)
  32. Wiersema SH, van Hooff HJA, Migchelsen LAA, Steultjens MPM. Reliability of the KT1000 arthrometer and the Lachman test in patients with an ACL rupture. *Knee.* 2008;15:107–110. doi:[10.1016/j.knee.2008.01.003](#)
  33. Hattori K, Ogawa M, Tanaka K, Matsuya A, Uematsu K, Tanaka Y. Can joint sound assess soft and hard endpoints on the Lachman test? *Biomed Mater Eng.* 2016;27:111–118. PubMed ID: [27175472](#)
  34. Cicchetti DV. Guidelines, criteria, and rules of thumb for evaluating normed and standardized assessment instruments in psychology. *Psychol Assess.* 1994;6(4):284–290. doi:[10.1037/1040-3590.6.4.284](#)